


PAJUNK®



InfiltraLong System

Continuous Wound Infiltration

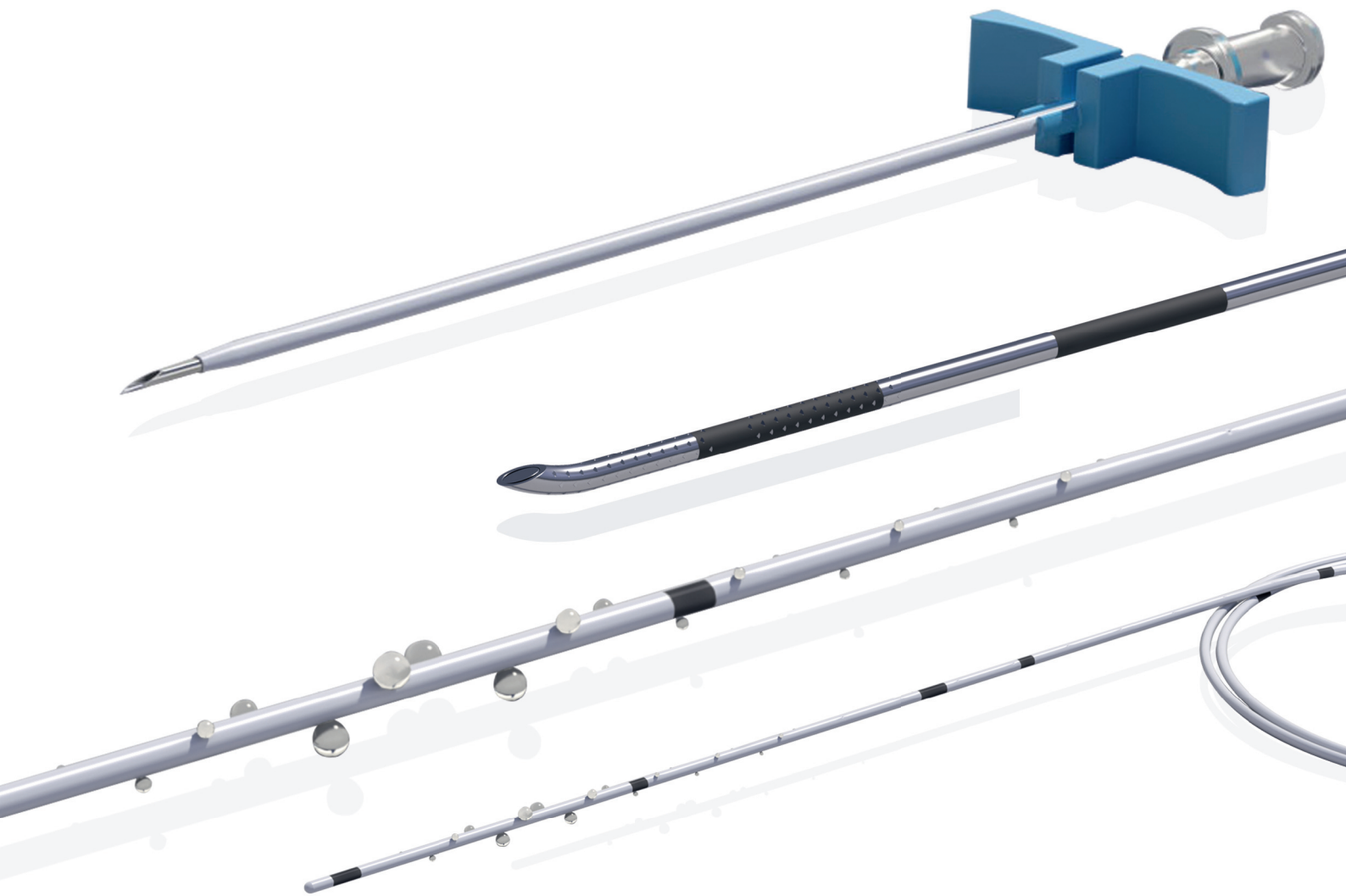

InfiltraLong System


The InfiltraLong catheter developed by Pajunk is a self-contained system for continuous infiltration analgesia at the surgical site. The design of the catheter allows a smooth and uniform targeted administration of local anaesthetic. This allows to optimise patient analgesia and improve outcomes while reducing opioid consumption.

The InfiltraLong catheter is the only wound infiltration catheter that can be combined with all commercial pumps. The wide range of catheter and perforation lengths available makes it suitable for many different requirements and areas of use. In comparison with three other commercial catheters the InfiltraLong performed extremely well because of its numerous fine infusion perforations, which are closely and evenly spaced in the perforation segment. This ensures an even distribution of the anaesthetic and a earlier mobilisation.



1. Catheter tip

- Transparent material
- Closed tip with lateral eyes

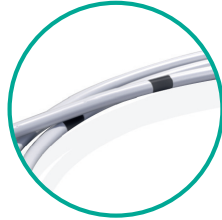

2. Depth Markings

- Graduation on the first 30cm
- Determines the exact position of catheter



InfiltraLong Set (with spilt cannula)



Item no.	Needle	InfiltraLong Catheter	Kit includes	PU
001158-00Z	Split cannula with puncture needle	19G x 420mm 15 openings along the first 40mm with helical coil	Filter FixoLong	10
001158-20B	Split cannula with puncture needle	19G x 600mm 60 openings along the first 150mm with helical coil	Filter FixoLong	10
001158-20D	Split cannula with puncture needle	19G x 600mm 45 openings along the first 100mm with helical coil	Filter FixoLong	10
001158-30C	Split cannula with puncture needle	19G x 600mm 88 openings along the first 220mm with helical coil	Filter FixoLong	10

InfiltraLong Set (with Tuohy needle)

001159-20D	Tuohy 17G x 90mm	19G x 600mm 15 openings along the first 40mm with helical coil	Clamp adaptor Filter & FixoLong	10
------------	------------------	---	------------------------------------	----

InfiltraLong Sono (Tuohy Sono needle)

001185-20B	Tuohy Sono 17G x 90mm	19G x 600mm 15 openings along the first 40mm with helical coil	Clamp adaptor Filter & FixoLong	10
------------	-----------------------	---	------------------------------------	----

Y-connector tube

Item no.	Description	PU
001151-37W	Injecting tube with y-connection for the simultaneous use of 2 catheters (Luer lock ends)	10



FixoCath

001161-37Z	Catheter fixation device for patient comfort	10
------------	--	----



PAJUNK® GmbH Medizintechnologie
Global Headquarters
Karl-Hall-Strasse 1
78187 Geisingen, Germany
Tel.: +49 7704 9291-0
Fax: +49 7704 9291-600
info@pajunk.com
pajunk.com

Do you have any questions? Our customer service & sales support will be glad to assist you.

Products distributed in Australia and New Zealand by

SURIMEX
www.surimex.com.au

15 Redland Drive
Mitcham VIC 3132
Australia

T (613) 8833 3804
E sales@surimex.com.au